

# DrawAlignR: An interactive tool for across run chromatogram alignment visualization

---

Justin Sing, Shubham Gupta, Arshia Mahmoodi, and Hannes Röst

# Contents

---

## 1. Data Preparation

- Performing Full OpenSwathWorkflow
  - Generating Library / Targeted Assay
  - Running an OpenSwathWorkflow
  - Scoring with PyProphet
- Required Input Data for DrawAlignR
- Downloading Tutorial Dataset

## 2. DrawAlignR User Manual

- Getting Started
- Performing Alignment
- Graphical Elements

# Data Preparation

---

DrawAlignR requires three types of files that are generated through an OpenSWATH Workflow using OpenMS tools.

For more information on OpenMS, please visit <https://www.openms.de/>

For more information on the OpenSWATH workflow, please visit <http://openswath.org/en/latest/>

Before running an OpenSwathWorkflow, you will need to ensure you have the correct input formats. You can convert all raw MS data to standard mzML/mzXML format using [ProteoWizard-MSConvert](#).

## Library / Assay Generation

### Generating the Target Assay Library

To generate the Targeted Assay library you need to run two OpenMS tools. [OpenSwathAssayGenerator](#) and [OpenSwathDecoyGenerator](#). The former will generate assays for each target peptide in the dda spectral library, the latter will append decoys to the targeted assay library from FDR estimation.

#### OpenSwathAssayGenerator

This tool has several parameters that you can adjust to tweak your target library. To understand what the parameters are, you can visit [https://abibuilder.informatik.uni-tuebingen.de/archive/openms/Documentation/release/latest/html/OpenMS\\_tutorial.html](https://abibuilder.informatik.uni-tuebingen.de/archive/openms/Documentation/release/latest/html/OpenMS_tutorial.html). This includes the documentation for the rest of OpenMS tools.

```
$ ./build/openms_build/bin/OpenSwathAssayGenerator -in
./dda_spectral_library/psgs_standard_consensus_Phospho_pep.TraML -out
./psgs_phospho_optimized.TraML
```

#### OpenSwathDecoyGenerator

To estimate true signal (true peptides at feature x), we need to differentiate/separate true positives from false positives. In order to do this, we use the decoy model that puts in false signal/peptides. This tool also has several different flags that you can set, such as different decoy methods: [shuffling](#), [reverse](#), [pseudo-reverse](#), etc. and switchKR which switches lysine/arginine ends. You can also adjust the ratio of decoy:targets.

```
$ ./build/openms_build/bin/OpenSwathDecoyGenerator -in
./psgs_phospho_optimized.TraML -out ./psgs_phospho_optimized_decoys.TraML
```

## TargetedFileConverter

Newer versions of OpenMS have moved towards sqlite data structures, due to more efficient data storage usage. You should convert the TraML target assay library generated in the previous step to the **pqp** sqlite format.

```
$ ./build/openms_build/bin/TargetedFileConverter -in  
psgs_phospho_optimized_decoys.TraML -out psgs_phospho_optimized_decoys.pqp
```

## OpenSwathWorkflow

Now that you have generated the targeted assay library, you are ready to run OpenSwathWorkflow.

### Default OpenSwathWorkflow

File names below are used as an example, make the output names more understandable.

**NOTE:** You need to run the OpenSwathWorkflow for each raw MS run file you have.

```
$ ./build/openms_build/bin/OpenSwathWorkflow \  
-in ./run1.mzXML \ ## Raw mzXML data file  
-tr ./lib/psgs/pqp \ ## Targeted assay library file  
-tr_irt ./IRT.TraML \ ## IRT peptides library file  
-out_osw ./run1_MSDATA_RESULTS.osw \ ## name of output file, output  
can be .tsv file or .osw (sqlite) file.  
-out_chrom run1.chrom.mzML \ ## if you want to generate  
chromatogram output file. Can be .mzML or sqMass (sqlite).  
-batchSize 1000 \ ## Size in which to split the data up into for  
extraction and scoring  
-readOptions cacheWorkingInMemory \ ## read the data on the fly, or  
read from cached in memory  
-tempDirectory ~/scratch/ \ ## temporary directroy to store cached  
files  
-threads 6 \ ## number of threads to use for computation  
-debug 10 ## debug level, debug 10 will spit out verbose  
information, helpful for debugging
```

## Pyprophet

After running OpenSwathWorkflow you want to score the results (extracted peak groups from OpenSwathWorkflow) and estimate the FDR distribution. To do that you can use [pyprophet](#).

Pyprophet contains several steps of it's own.

### Pyprophet Merge

You first want to merge each individual osw run file into one osw file, this make things easier for running the workflow and getting global estimates. You can also run it on each file individually, but generally we run the workflow on a merged.osw file.

```
pyprophet merge --template=./lib/psgs_phospho_optimized_decoys.pqp --  
out=merged.osw $(find ./results/ -maxdepth 2 -type f -name  
*_MSDATA_RESULTS.osw)
```

The command above will call the pyprophet merge function, and merge all the osw files (that are found using the bash find function to return all the files matching "\*\_MSDATA\_RESULTS.osw"), into one osw file using the library file as a template.

### Pyprophet Scoring

After you have a merged osw file, you want to perform the semi-supervised scoring algorithm.

Generally you would only need to score on the MS2 level. However, you can also use an integrated scoring on the --level=**ms1ms2**, this will integrate the scores between MS1 features and MS2 features.

```
pyprophet score --in=./merged.osw --level=ms2 --ss_num_iter=10 --  
xeval_num_iter=10 --threads=10
```

The command above will score on the merged\_runs.osw file, using the MS2 features, using the semi-supervised learning on 10 iterations (this is the default) and performing 10 cross-validations (this is the default) using 10 threads.

## Required Input Data for DrawAlignR

### Directory Structure

To make things easier, and to have a contained experiment, DrawAlignR utilizes a structured directory as input.

The top level experiment directory should contain sub-directories:

- containing the merged pyprophet scored **osw** file obtained from the OpenSwathWorkflow and pyprophet scoring workflow,
- the extracted chromatogram files (**mzML** or **sqMass**) obtained from OpenSwathWorkflow, and
- the **pqp** library assay file obtained from OpenSwathAssayGenerator and OpenSwathDecoyGenerator.

DrawAlignR, can alternatively also allow for individual file input without having a structured working directory, however, it is suggested to have one.

Please see below for an example structured working directory.

**An example Working Directory should be ideally structured and named as below:**

```
/Project_Working_Directory
|__ /mzml
|   |__ /run0.chrom.mzML
|   |__ /run1.chrom.mzML
|__ /osw
|   |__ /merged.osw
|__ /pqp
|   |__ /assay_library.pqp
|__ /sqmass(Optional, either mzml or sqmass format is acceptable)
|   |__ /run1.chrom.sqMass
|   |__ /run1.chrom.sqMass
```

### Data Size

Depending on how complex your experiments are and how many runs you may have, the extracted chromatogram files may be large files. The larger your files are, the longer it may take for DrawAlignR to cache all chromatographic data into memory, this is something to keep in mind.

## Downloading Tutorial Dataset

We have example datasets hosted on PeptideAtlas [PASS01520](#)

This dataset contains two datasets, which both have a full dataset and a smaller data subset:

- Spyogenes
  - Spyogenes\_Full
  - Spyogenes\_Small\_Subset
- Synthetic\_Phosphopeptide\_Dataset
  - Synthetic\_Dilution\_Phosphoproteomics\_Full
  - Synthetic\_Dilution\_Phosphoproteomics\_Small\_Subset

**NOTE:** Using the smaller data subsets, will be quicker for tests, and alignment parameter checks, since loading full extracted chromatogram data may take longer to cache into memory.

The easiest and most straight forward way of downloading the data is to use your browser's FTP mode:  
`ftp://PASS01520:HE7445u@ftp.peptideatlas.org/`

### Retrieving Data on MacOS

MacOS is capable of establishing a connection to a remote FTP server using the native Finder application.

- Open a **Finder Window**
- From the main menu, select **Go -> connect to Server**
- In the pop-up dialog box, type `ftp://PASS01520:HE7445u@ftp.peptideatlas.org/` in the *Server Address:* field
- Press **Connect**
- enter Username/Password
  - Username: PASS01520
  - Password: HE7445u

Alternatively, you can type the ftp address into a safari browser.

### Retrieve Data on Windows

On Windows, you can use the Windows file manager/ File Explorer to connect to a remote FTP server.

- Open a **File Explorer**
- Right click on **This PC / Computer** \*\* Select **Add a network location**
- Go through pop-up wizard and select **Choose a custom network location**
- In the "Specify the location of your website" dialog, enter  
`ftp://PASS01520:HE7445u@ftp.peptideatlas.org/`
- enter Username/Password
  - Username: PASS01520
  - Password: HE7445u

## Retrieve Data on Linux/ Ubuntu

On Ubuntu, you can use nautilus file manager to connect to a remote FTP server.

- Click on **Other Locations**
- In the bottom field to the right of *Connect to Server\**, enter  
`ftp://PASS01520:HE7445u@ftp.peptideatlas.org/`
- enter Username/Password
  - Username: PASS01520
  - Password: HE7445u



# Using/ Calling DrawAlignR

---

The visualization tool can be called via two methods, either via installing the repository using devtools github package installer or by calling shiny::runGitHub.

## Calling tool via Installing DrawAlignR

```
### Install devtools if not already installed.  
## install.packages("devtools")  
## Install DrawAlignR from Roestlab Github Repository  
devtools::install_github("Roestlab/DrawAlignR")  
## To run the tool  
DrawAlignR::runDrawAlignR()
```

OR

## Calling tool via Github repo

```
### Install shiny  
## install.packages("shiny")  
shiny::runGitHub(repo = "Roestlab/DrawAlignR/", username = "Roestlab",  
  subdir = "inst/shiny-script")
```

# User Interface

---

There are three major tabs in the left side pannel:

- **General Settings**
  - These are general settings for uploading chromatogram files (.mzML or .sqMass), library assay files (.ppq), OpenSwathWorkflow results file (.osw). The user can also set the working directory that contains sub-directories for osw and mzml files.
  - The user selects the peptide and charge state to visualize.
  - The user selects how many plots to show for each chromatogram run file supplied.
  - The user can select the alignment option to perform an alignment for the selected peptide.
  - The user can visualize the reference plot, experiment plot and the experiment aligned plot.
- **Alignment Settings**
  - The user can change various alignment parameters
- **Plot Settings**
  - The user can change various plot visualization settings

/media/justincsing/ExtraDrive1/Documents2/Roest\_Lab/Github/DrawAlignR/inst/shiny-script - Shiny

http://127.0.0.1:4555 | Open in Browser | Publish

**DrawAlignR** Ver: 0.1.0

**Error: A Chromatogram file was not supplied or not found**

General Settings | Alignment Settings | Plot Settings

Working Directory ☒ ?

Set Working Directory

Set Working Directory

Peptide Name

Peptide Charge

☐ Plot Aligned

Select Reference Run for Alignment

Experiment(s) to Align

☐ Show Original Peak Annotation

Toggle to switch between using a working directory or supplying individual files

Set Working Directory  
- osw sub-directroy  
- mzml sub-directory etc.

User click button  
or  
User input text box

Peptide Selection

Peptide Charge Selection

Perform Alignment

Run to use as Reference

Run to use as Experiment

Option to show original peak annotation

User Searchable Drop-down list

User Searchable Drop-down list

User Searchable Drop-down list

User Searchable Drop-down list

User Searchable Drop-down list

User Searchable Drop-down list

User Checkbox

User Checkbox


# Tutorial for Performing Alignment

## Set Working Directory

Use the Set Working Directory button to set the working directory that contains an mzml folder with .chrom.mzml files and an osw file with a merged.osw file, and optionally a library file in the .ppq format. Or you can directly enter the path to the working directory using the input textbox area. DrawAlignR will parse the working directory and look for `../osw`, `../mzml` ( or `../sqmass` ), and `../ppq` subfolders containing their corresponding file types. You can alternatively toggle the **Working Directory** button to enter chromatogram, library and osw files separately. **Note:** For alignment it is preferred if you enter a working directory.

An example Working Directory should be ideally structured and named as below:

```
/Project_Working_Directory
|__/_mzml
|   |__/_run0.chrom.mzML
|   |__/_run1.chrom.mzML
|__/_osw
|   |__/_merged.osw
|__/_ppq
|   |__/_assay_library.pqp
|__/_sqmass(Optional, either mzml or sqmass format is acceptable)
|   |__/_run1.chrom.sqMass
|   |__/_run1.chrom.sqMass
```

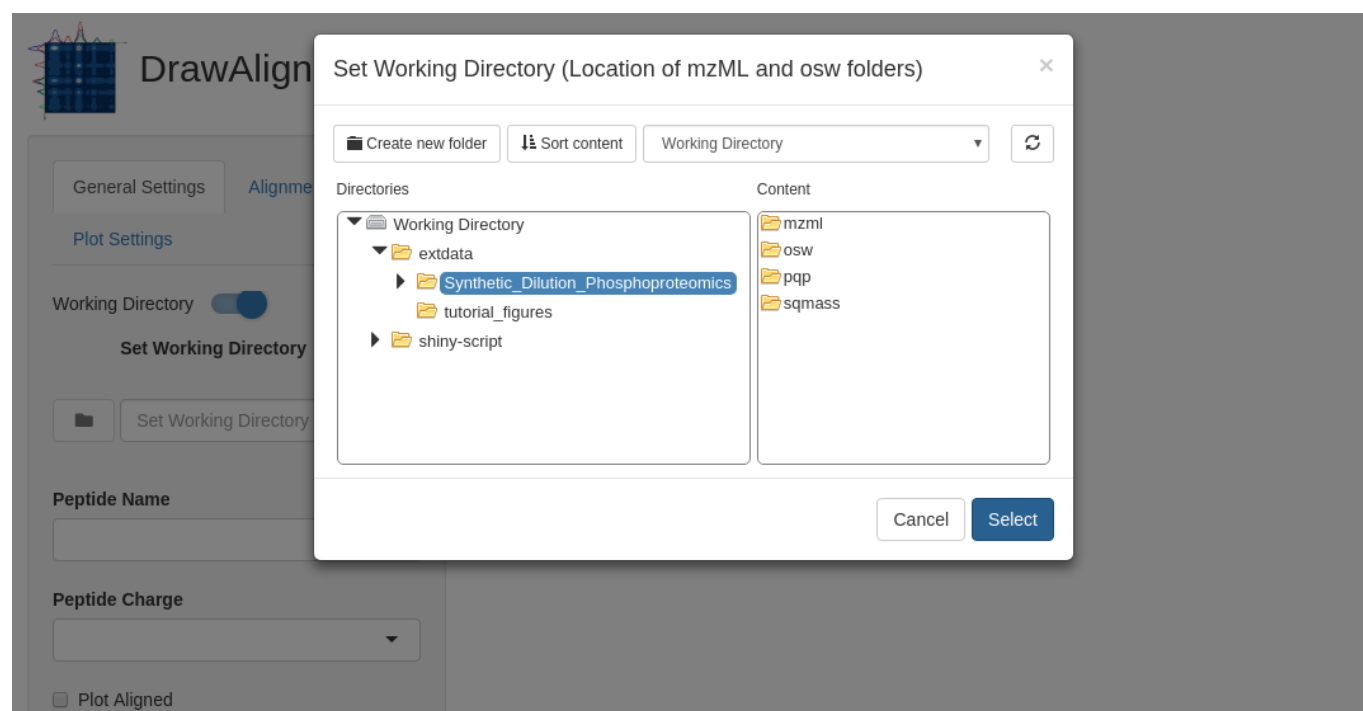

## Select a Peptide to visualize Chromatogram alignment

Choose which peptide you want to visualize using the Peptide dropdown list. The dropdown list is searchable, so you can easily search for a specific peptide to visualize. The list of peptides is extracted from either the input library file if available, or an osw file if available.

The screenshot displays the DrawAlignR web application interface. The browser address bar shows the URL `http://127.0.0.1:5681`. The application title is `/media/justincsing/ExtraDrive1/Documents2/Roest_Lab/Github/DrawAlignR/inst/shiny-script - Shiny`. The application logo is **DrawAlignR** Ver: 0.1.0.

**General Settings**

- Alignment Settings** (selected)
- Plot Settings**
- Working Directory**: `/media/justincsing/ExtraDrive1/Do`
- Set Working Directory** button
- ☒ **Use mzml** ☐ **Use sqmass**
- Peptide Name** dropdown menu with a search bar and a list of peptides including: ANS(UniMod:21)SPTTNIDHLK(UniMod:2), ANSS(UniMod:21)PTTNIDHLK(UniMod:2), ANSSPT(UniMod:21)TNIDHLK(UniMod:2), ES(UniMod:21)TAEPDSLRS(UniMod:267), EST(UniMod:21)AEPDSLRS(UniMod:267), ES(UniMod:21)T(UniMod:21)AEPDSLRS, ESTAEPDS(UniMod:21)LSR(UniMod:267), KDS(UniMod:21)NTNVLK(UniMod:250).
- ☐ **Plot Aligned**
- Select Reference Run for Alignment**: `chludwig_K150309_001b_SW_1_64`
- Experiment(s) to Align** list: chludwig\_K150309\_002b\_SW\_1\_32, chludwig\_K150309\_003b\_SW\_1\_24, chludwig\_K150309\_004b\_SW\_1\_16, chludwig\_K150309\_005b\_SW\_1\_12, chludwig\_K150309\_006b\_SW\_1\_8, chludwig\_K150309\_007b\_SW\_1\_6, chludwig\_K150309\_008\_SW\_1\_4, chludwig\_K150309\_009\_SW\_1\_3, chludwig\_K150309\_010\_SW\_1\_2, chludwig\_K150309\_011\_SW\_1\_1point5, chludwig\_K150309\_012\_SW\_1\_1, chludwig\_K150309\_013\_SW\_0.
- ☐ **Show Original Peak Annotation**

**Chromatogram Plots**

Three empty chromatogram plots are displayed, each with a title bar showing `Run: 001b_1_64 | Precursor: | Peptide: | Charge: | m/z:`. The plots have x-axes from 0 to 1 and y-axes from 0 to 1. Each plot has a toolbar with icons for zooming, panning, and other plot controls.

# Select a Charge State for Selected Peptide

Choose which charge state to visualize for the selcted peptide.

General Settings

Alignment Settings

Plot Settings

Working Directory

Set Working Directory

/media/justincsing/ExtraDrive1/Do

Use mzml

Use sqmass

Peptide Name

EST(UniMod:21)AEPDSLRSR(UniMod:267)

Peptide Charge

2

2:3

2

3

Run 9

Run 10

Run 11

Run 12

Run 13

Plot Aligned

Select Reference Run for Alignment

chludwig\_K150309\_001b\_SW\_1\_64

Experiment(s) to Align

chludwig\_K150309\_002b\_SW\_1\_32

chludwig\_K150309\_003b\_SW\_1\_24

chludwig\_K150309\_004b\_SW\_1\_16

chludwig\_K150309\_005b\_SW\_1\_12

chludwig\_K150309\_006b\_SW\_1\_8

chludwig\_K150309\_007b\_SW\_1\_6

chludwig\_K150309\_008\_SW\_1\_4

chludwig\_K150309\_009\_SW\_1\_3

chludwig\_K150309\_010\_SW\_1\_2

chludwig\_K150309\_011\_SW\_1\_1point5

chludwig\_K150309\_012\_SW\_1\_1

chludwig\_K150309\_013\_SW\_0

Show Original Peak Annotation

EST(UniMod:21)AEPDSLRSR(UniMod:267)

Run: 001b\_1\_64 | m/z: 641.2653

Int

1200

1000

800

600

400

200

0

RT

1000

1200

1400

Detecting

1068\_1+\_y6\_684.3551

1069\_1+\_y4\_472.2753

1070\_1+\_y7\_813.3976

1071\_1+\_y8\_884.4348

1072\_1+\_y5\_587.3023

1073\_1+\_y1\_185.1272

EST(UniMod:21)AEPDSLRSR(UniMod:267)

Run: 002b\_1\_32 | m/z: 641.2653

Int

1200

1000

800

600

400

200

0

RT

1000

1200

1400

Detecting

1068\_1+\_y6\_684.3551

1069\_1+\_y4\_472.2753

1070\_1+\_y7\_813.3976

1071\_1+\_y8\_884.4348

1072\_1+\_y5\_587.3023

1073\_1+\_y1\_185.1272

EST(UniMod:21)AEPDSLRSR(UniMod:267)

Run: 003b\_1\_24 | m/z: 641.2653

Int

1200

1000

800

600

400

200

0

RT

1000

1200

1400

Detecting

1068\_1+\_y6\_684.3551

1069\_1+\_y4\_472.2753

1070\_1+\_y7\_813.3976

1071\_1+\_y8\_884.4348

1072\_1+\_y5\_587.3023

1073\_1+\_y1\_185.1272

# Select Reference Run

Choose which chromatogram file to use as the reference run. These are set through the searchable dropdown lists, which extracts the filenames from the supplied chromatogram files without the .chrom.mzml extension

General Settings

Alignment Settings

Plot Settings

Working Directory

Set Working Directory

/media/justincsing/ExtraDrive1/Do

☒ Use mzml

☐ Use sqmass

Peptide Name

EST(UniMod:21)AEPDSLRSR(UniMod:267)

Peptide Charge

2

Choose Runs to Display:

☒ Run 1

☒ Run 2

☒ Run 3

☒ Run 4

☒ Run 5

☒ Run 6

☒ Run 7

☒ Run 8

☒ Run 9

☒ Run 10

☒ Run 11

☒ Run 12

☒ Run 13

☐ Plot Aligned

Select Reference Run for Alignment

ch

chludwig\_K150309\_013\_SW\_0

chludwig\_K150309\_008\_SW\_1\_4

chludwig\_K150309\_009\_SW\_1\_3

chludwig\_K150309\_010\_SW\_1\_2

chludwig\_K150309\_012\_SW\_1\_1

chludwig\_K150309\_006b\_SW\_1\_8

chludwig\_K150309\_007b\_SW\_1\_6

chludwig\_K150309\_001b\_SW\_1\_64

chludwig\_K150309\_007b\_SW\_1\_6

chludwig\_K150309\_008\_SW\_1\_4

chludwig\_K150309\_009\_SW\_1\_3

chludwig\_K150309\_010\_SW\_1\_2

chludwig\_K150309\_011\_SW\_1\_1point5

chludwig\_K150309\_012\_SW\_1\_1

chludwig\_K150309\_013\_SW\_0

☐ Show Original Peak Annotation

EST(UniMod:21)AEPDSLRSR(UniMod:267)

Run: 001b\_1\_64 | m/z: 641.2653

Detecting

1068\_1+\_y6\_684.3551

1069\_1+\_y4\_472.2753

1070\_1+\_y7\_813.3976

1071\_1+\_y8\_884.4348

1072\_1+\_y5\_587.3023

1073\_1+\_y1\_185.1272

EST(UniMod:21)AEPDSLRSR(UniMod:267)

Run: 002b\_1\_32 | m/z: 641.2653

Detecting

1068\_1+\_y6\_684.3551

1069\_1+\_y4\_472.2753

1070\_1+\_y7\_813.3976

1071\_1+\_y8\_884.4348

1072\_1+\_y5\_587.3023

1073\_1+\_y1\_185.1272

EST(UniMod:21)AEPDSLRSR(UniMod:267)

Run: 003b\_1\_24 | m/z: 641.2653

Detecting

1068\_1+\_y6\_684.3551

1069\_1+\_y4\_472.2753

1070\_1+\_y7\_813.3976

1071\_1+\_y8\_884.4348

1072\_1+\_y5\_587.3023

1073\_1+\_y1\_185.1272

## Select Experiment Run(s)

Choose which chromatogram file to use as the experiment run. These are set through the searchable dropdown lists, which extracts the filenames from the supplied chromatogram files without the .chrom.mzml extension

Shiny interface for DrawAlignR (Ver: 0.1.0) showing settings and three chromatogram plots.

**General Settings**

- Working Directory: ☐ ☒ Set Working Directory:
- ☒ Use mzml ☐ Use sqmass
- Peptide Name:
- Peptide Charge:
- Choose Runs to Display:
  - ☒ Run 1 ☒ Run 2 ☒ Run 3 ☒ Run 4
  - ☒ Run 5 ☒ Run 6 ☒ Run 7 ☒ Run 8
  - ☒ Run 9 ☒ Run 10 ☒ Run 11
  - ☒ Run 12 ☒ Run 13
- ☐ Plot Aligned
- Select Reference Run for Alignment:
- Experiment(s) to Align:
  - chludwig\_K150309\_001b\_SW\_1\_64
  - chludwig\_K150309\_002b\_SW\_1\_32
  - chludwig\_K150309\_003b\_SW\_1\_24
  - chludwig\_K150309\_004b\_SW\_1\_16
  - chludwig\_K150309\_005b\_SW\_1\_12
  - chludwig\_K150309\_006b\_SW\_1\_8
  - chludwig\_K150309\_007b\_SW\_1\_6
  - chludwig\_K150309\_008\_SW\_1\_4
  - chludwig\_K150309\_009\_SW\_1\_3
  - chludwig\_K150309\_010\_SW\_1\_2
  - chludwig\_K150309\_011\_SW\_1\_1point5
  - chludwig\_K150309\_012\_SW\_1\_1 |
- ☐ Show Original Peak Annotation

**Chromatogram Plots (RT vs Int):**

EST(UniMod:21)AEPDSLRSR(UniMod:267) Run: 001b\_1\_64 | m/z: 641.2653

EST(UniMod:21)AEPDSLRSR(UniMod:267) Run: 002b\_1\_32 | m/z: 641.2653


EST(UniMod:21)AEPDSLRSR(UniMod:267) Run: 003b\_1\_24 | m/z: 641.2653

Detecting

- 1068\_1+\_y6\_684.3551
- 1069\_1+\_y4\_472.2753
- 1070\_1+\_y7\_813.3976
- 1071\_1+\_y8\_884.4348
- 1072\_1+\_y5\_587.3023
- 1073\_1+\_y1\_185.1272

## Select Which Chromatogram Runs to Display

Choose which chromatogram runs to display by checking or unchecking the Run **n** checkbox.

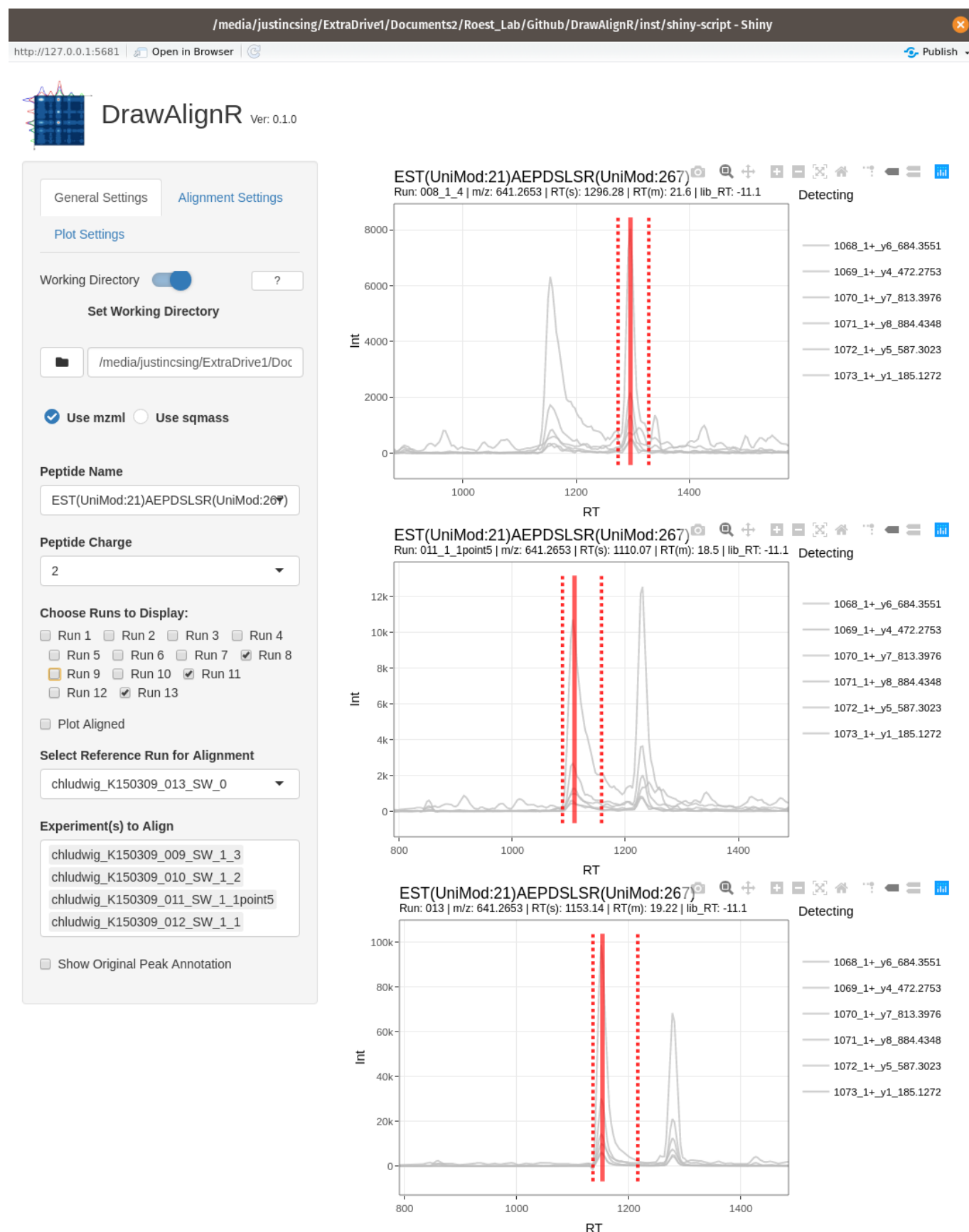




## Select Plot Align checkbox to perform alignment

Check the Plot Aligned checkbox to perform the alignment of the two runs for the selected peptide You can also use the

- Show Original Peak Annotation to show where the original peak was annotated before alignment.

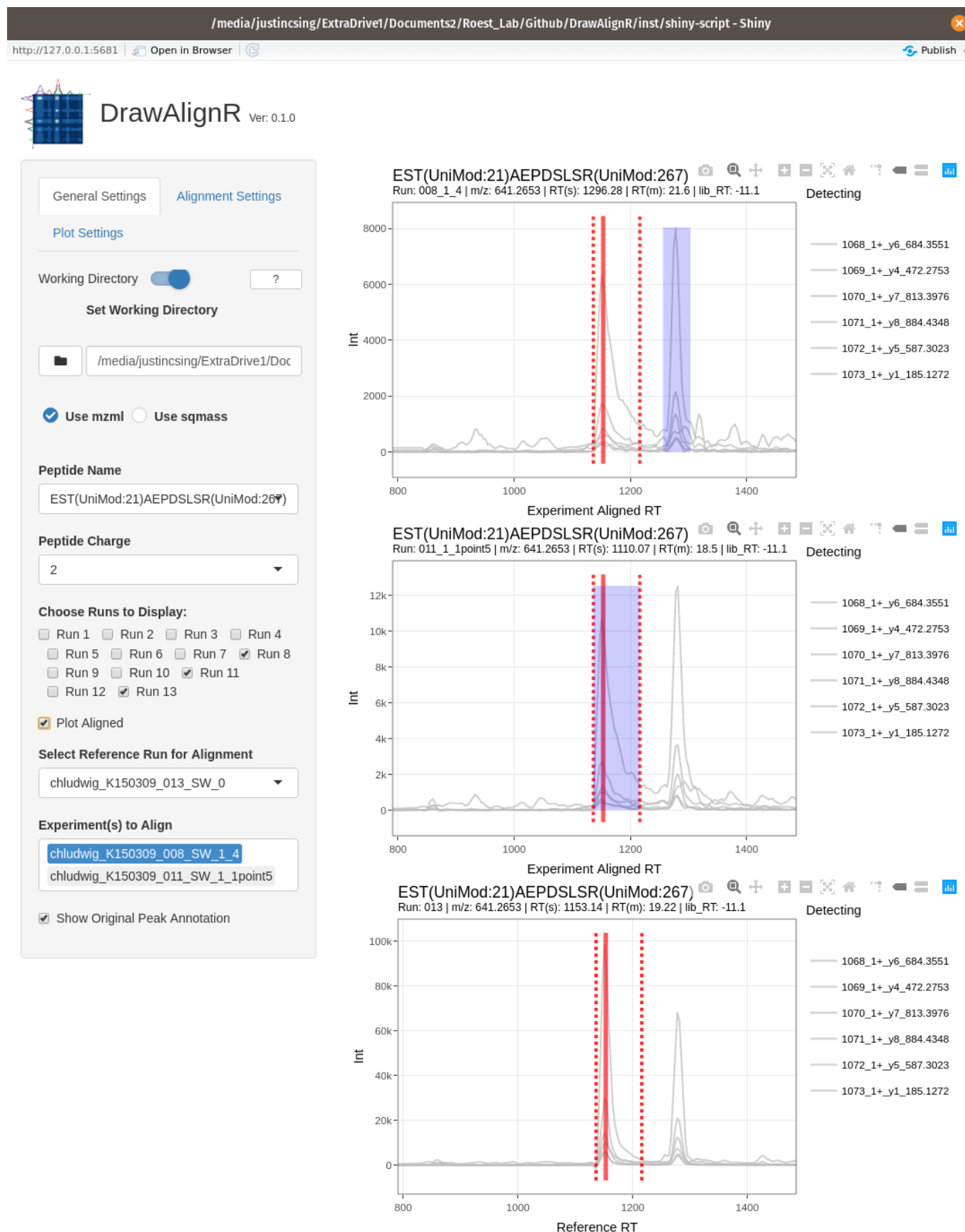


# Change Alignment Parameters

Select the Alignment settings tab to change various alignment settings




## Zoom into each chromatogram for further inspection

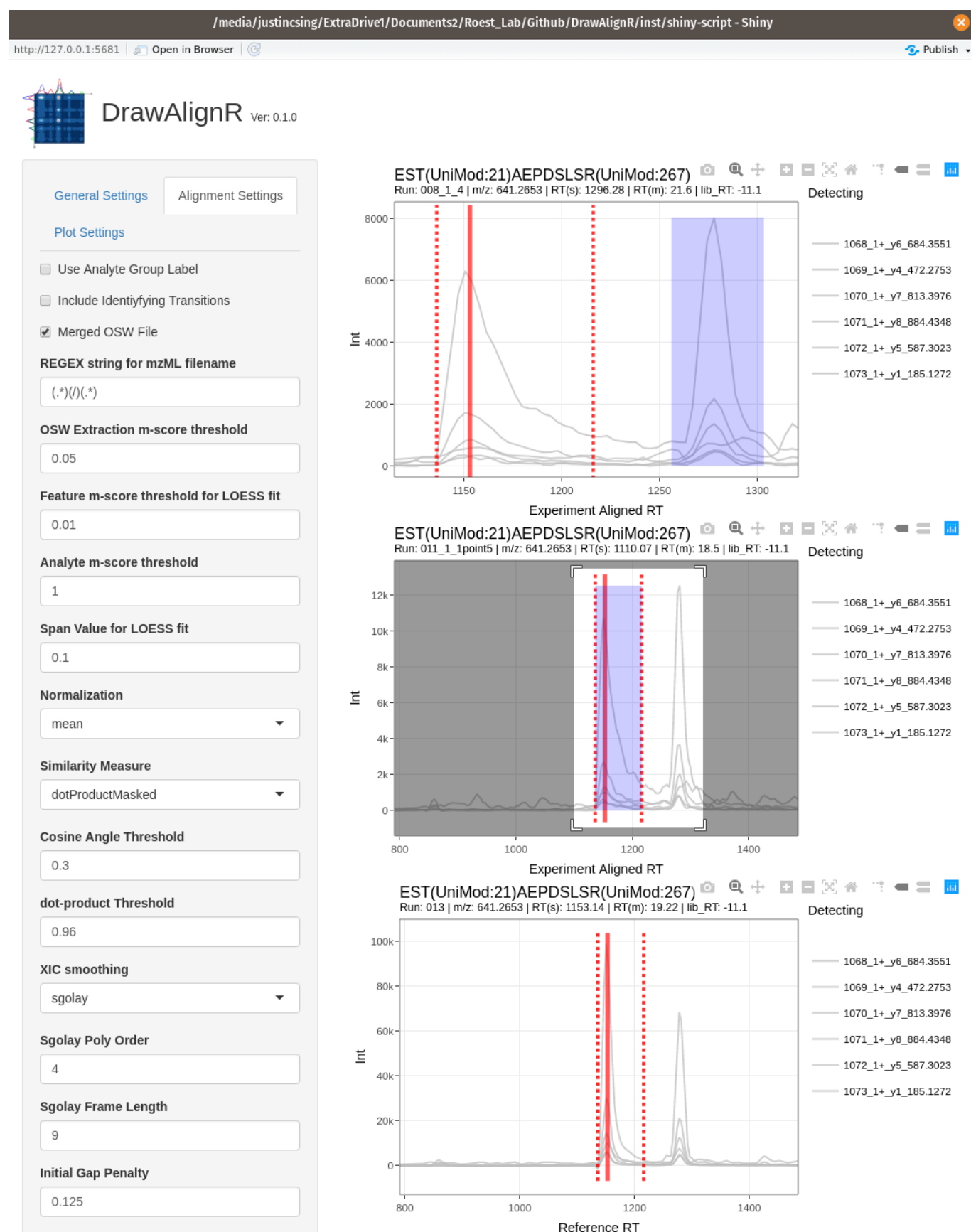


# Use the Hover-tooltip

You can hover over the chromatogram traces to see information such as Retention time and Intensity

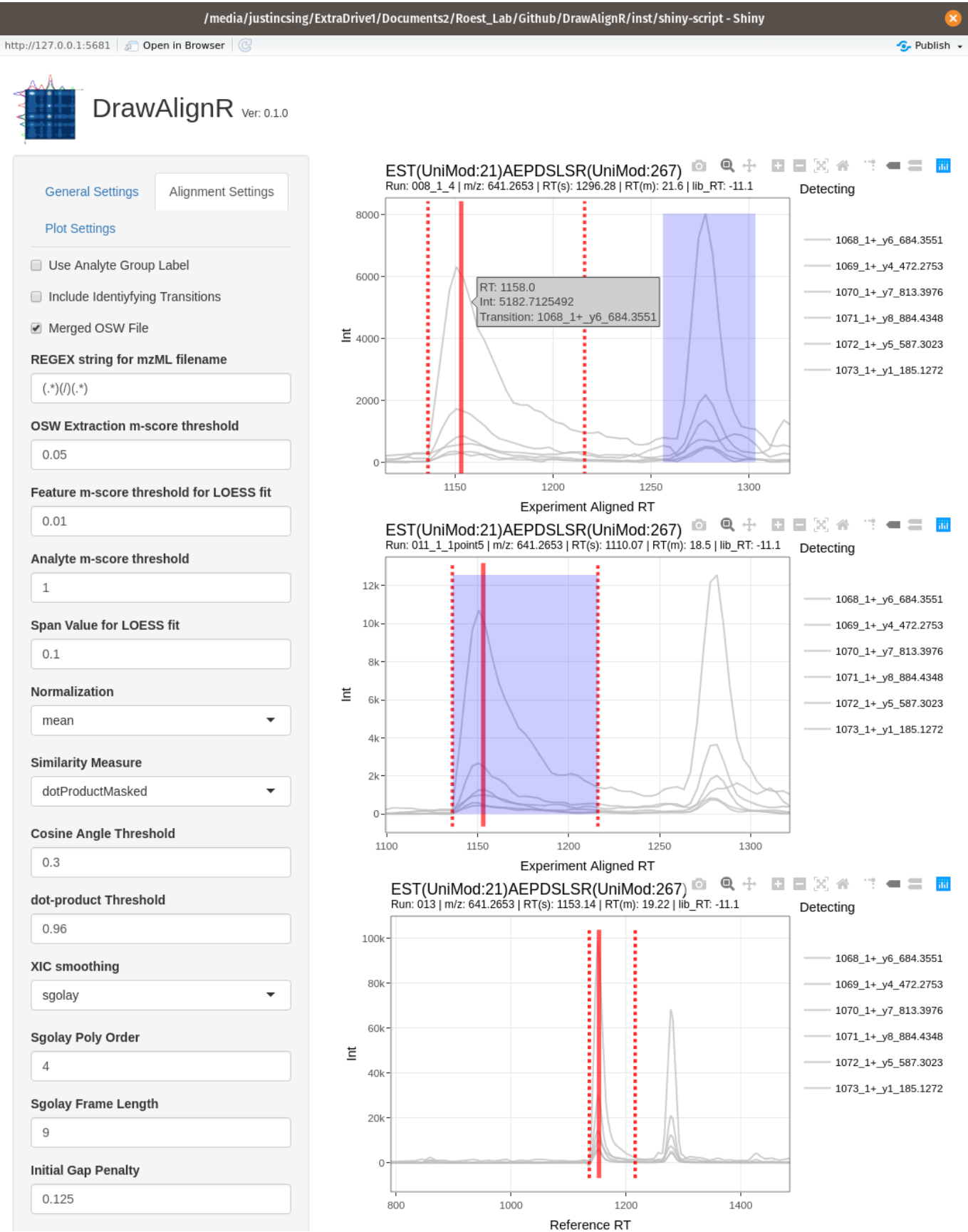

# Visualizing Chromatograms

The user can also just visualize the chromatograms alone without performing alignment to visually inspect each trace. If the user has an IPF dataset, they can visualize site-determining ions (unique identifying transitions) of the modified peptide.

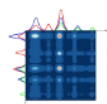

DrawAlignR Ver: 0.1.0

General Settings

Alignment Settings

Plot Settings

Working Directory

Chromatogram Path(s)

X

PQP Path

X

OSW Path

X

Peptide Name

ANS(UniMod:21)SPTTNIDHLK(UniMod:259)

Peptide Charge

2

Choose Runs to Display:

☒ Run 1 ☒ Run 2 ☒ Run 3

☐ Plot Aligned




Select Reference Run for Alignment

chludwig\_K150309\_013\_SW\_0

Experiment(s) to Align


chludwig\_K150309\_007b\_SW\_1\_6  
chludwig\_K150309\_004b\_SW\_1\_16

☐ Show Original Peak Annotation

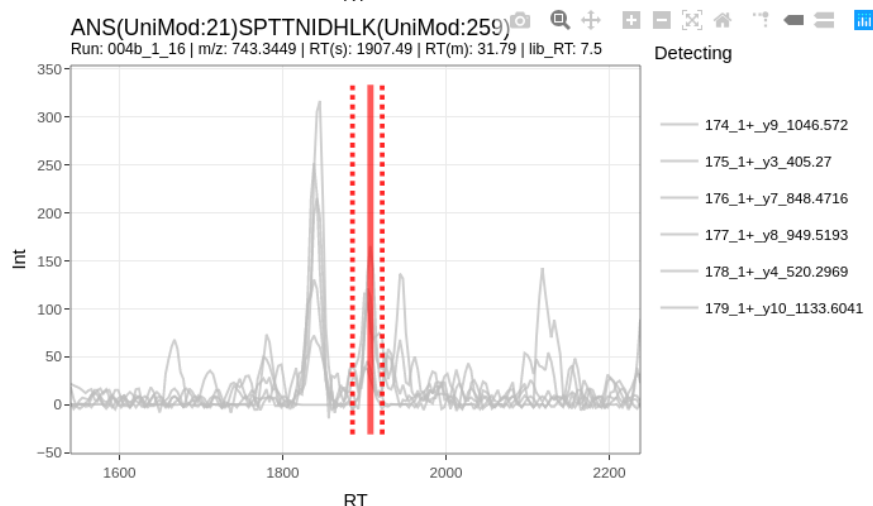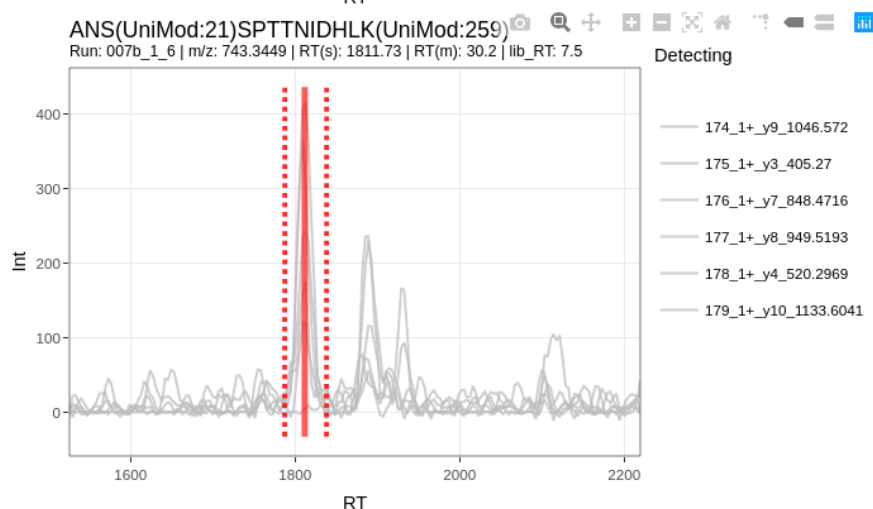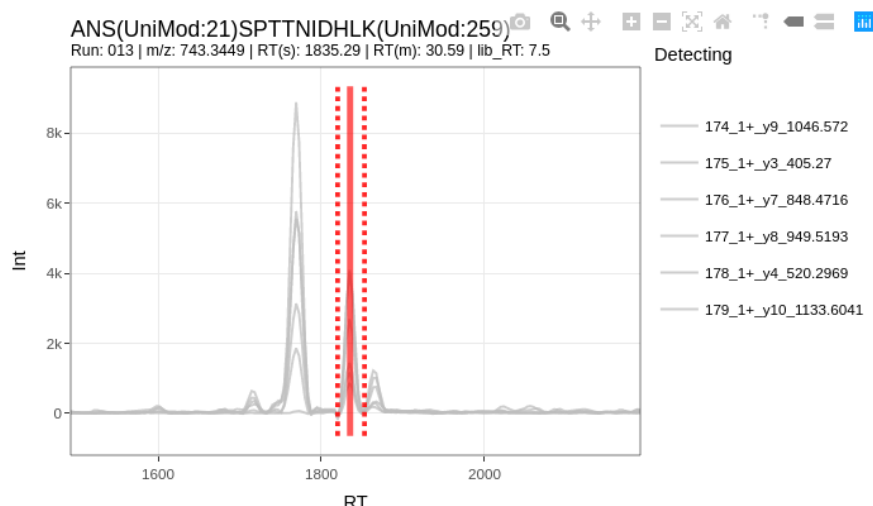

## Precursor, Detecting, Identifying

The user can choose to display the precursor trace, or the 6 detecting traces, or the unique identifying traces. The precursor trace is displayed in **black**, the detecting traces are displayed in a **light gray** and the unique identifying transitions are **colored**

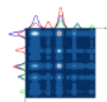

DrawAlignR Ver: 0.1.0

General Settings
Alignment Settings

Plot Settings




☒ Plot Precursor Trace  
☒ Plot Detecting Traces  
☒ Plot Unique Identifying Traces

Show n Identifying Traces

identifying y-ions


identifying b-ions

☒ Show Transition Scores (hover tooltip)  
☐ Show All Peak-Groups

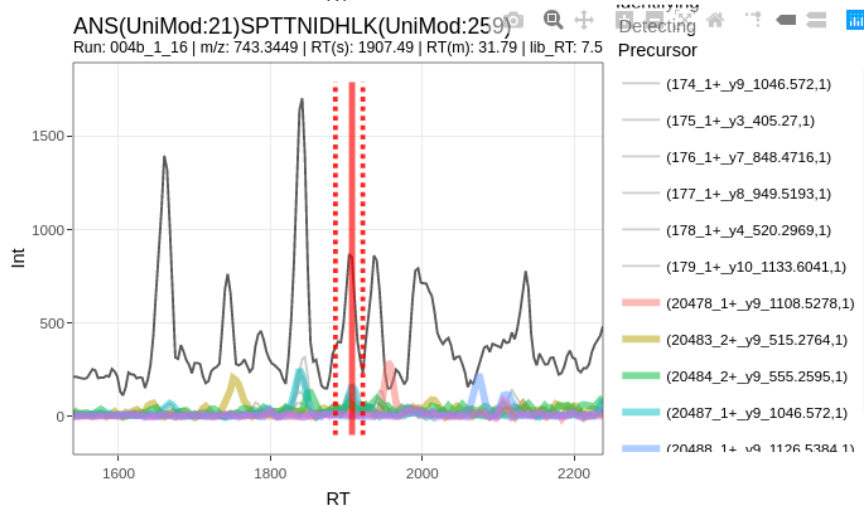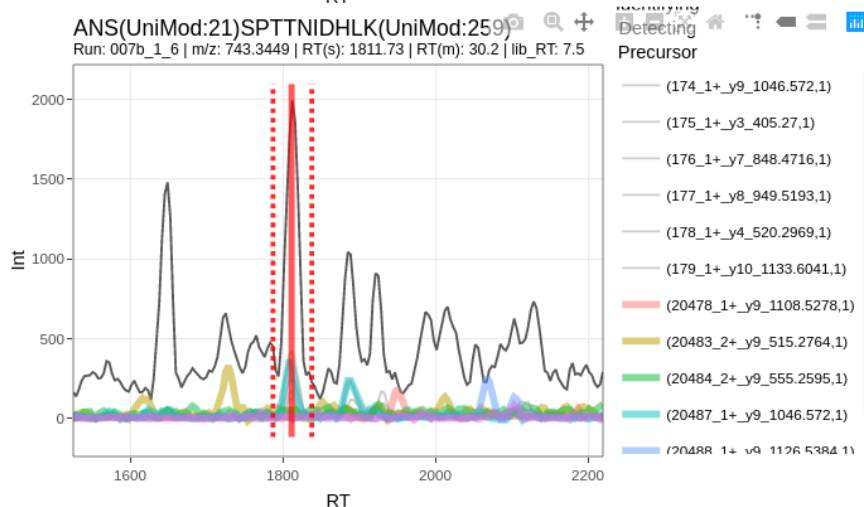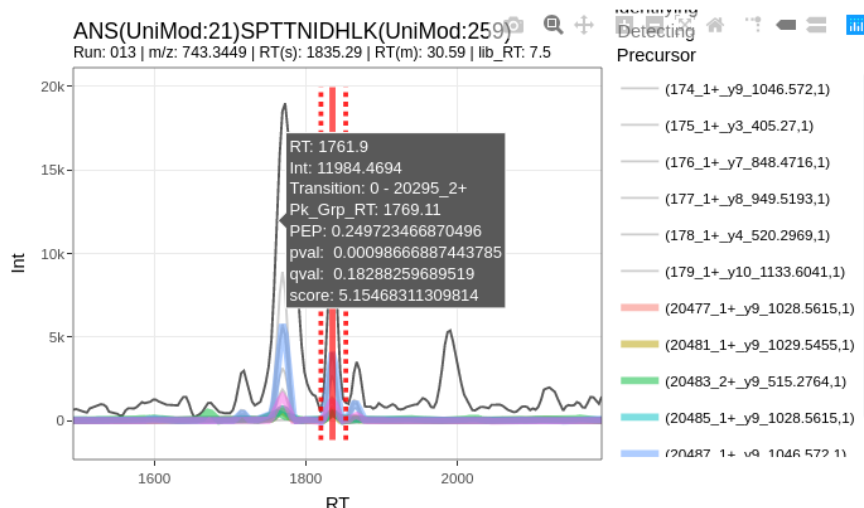

## Displaying Transition Scores

The user can hover of the traces to display the transition scores such as the transitions posterior error probability, q-value and score.

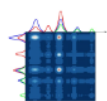

DrawAlignR Ver: 0.1.0

General Settings
Alignment Settings

Plot Settings


☒ Plot Precursor Trace  
☒ Plot Detecting Traces  
☒ Plot Unique Identifying Traces

Show n Identifying Traces

identifying y-ions

identifying b-ions


☒ Show Transition Scores (hover tooltip)  
☐ Show All Peak-Groups

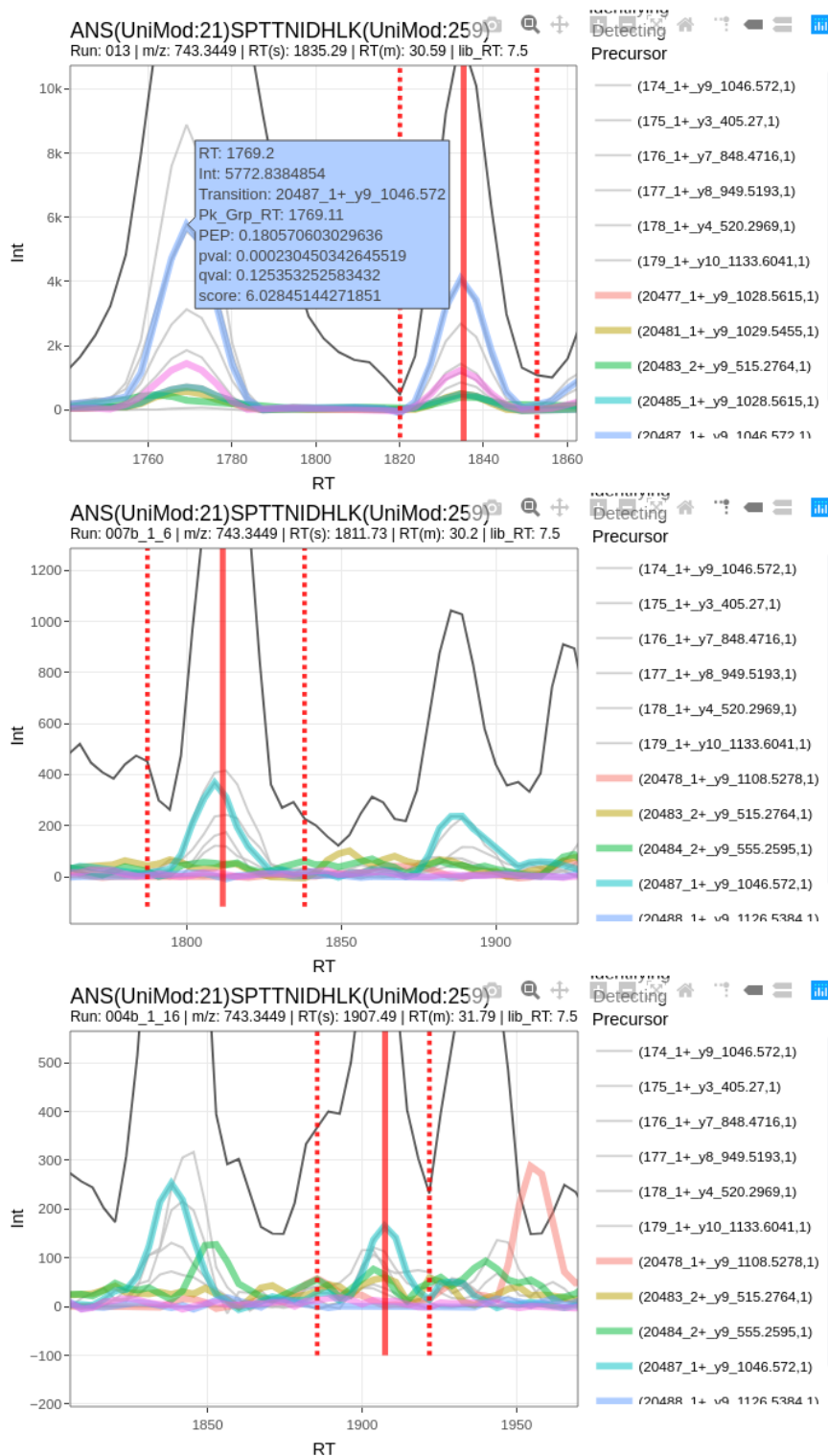



# Displaying all Peak-Group Ranks

The user can choose to display the other potential peak-group ranks found by OpenSWATH

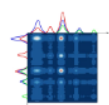

DrawAlignR Ver: 0.1.0

General Settings
Alignment Settings

Plot Settings

☒ Plot Precursor Trace

☒ Plot Detecting Traces

☒ Plot Unique Identifying Traces

Show n Identifying Traces

6

identifying y-ions




9,10

identifying b-ions

3

☒ Show Transition Scores (hover tooltip)


☒ Show All Peak-Groups

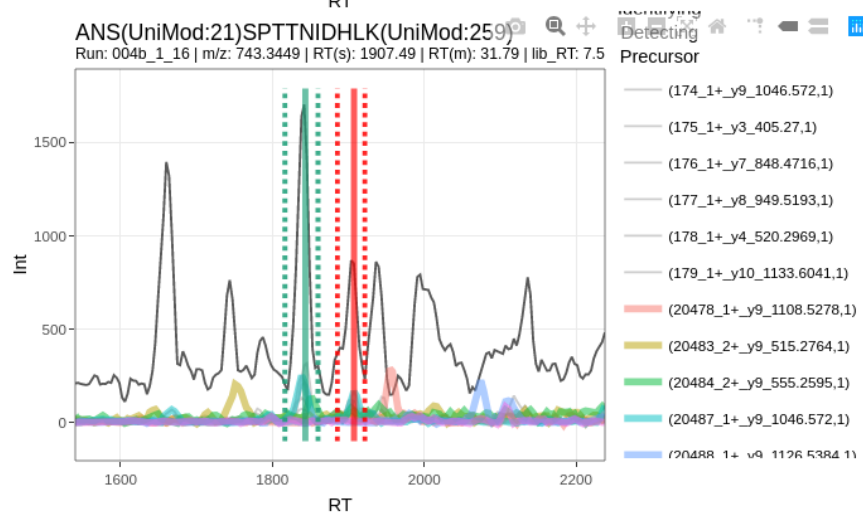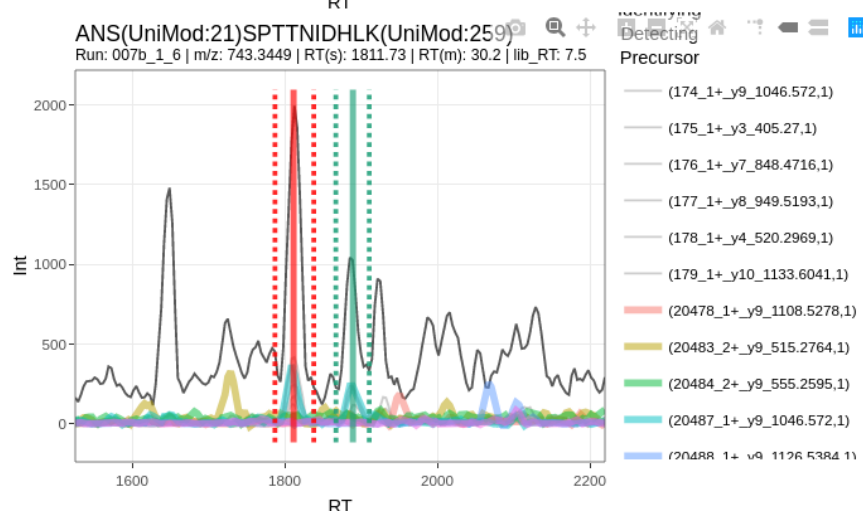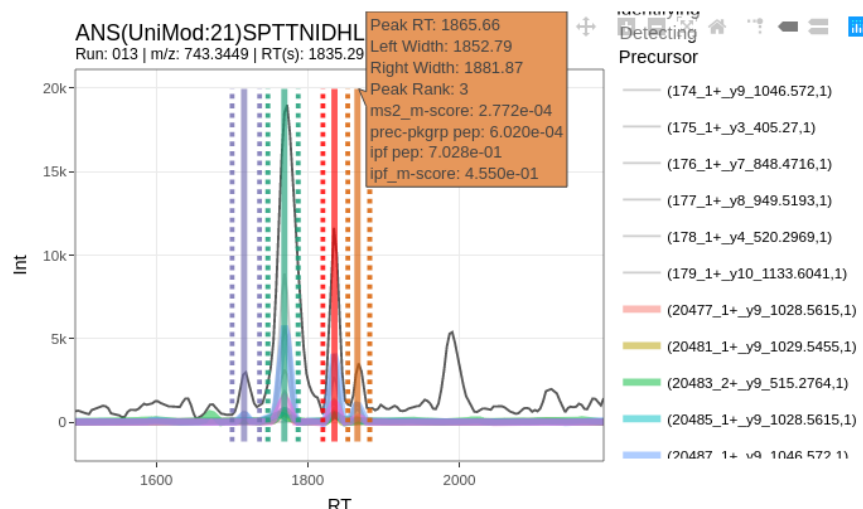



# Unselecting a Few Transitions to Display

The user can click on the legend to hide transitions they don't want to display

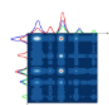

DrawAlignR Ver: 0.1.0

General Settings
Alignment Settings

Plot Settings




☒ Plot Precursor Trace  
☒ Plot Detecting Traces  
☒ Plot Unique Identifying Traces

Show n Identifying Traces

identifying y-ions


identifying b-ions

☒ Show Transition Scores (hover tooltip)  
☒ Show All Peak-Groups

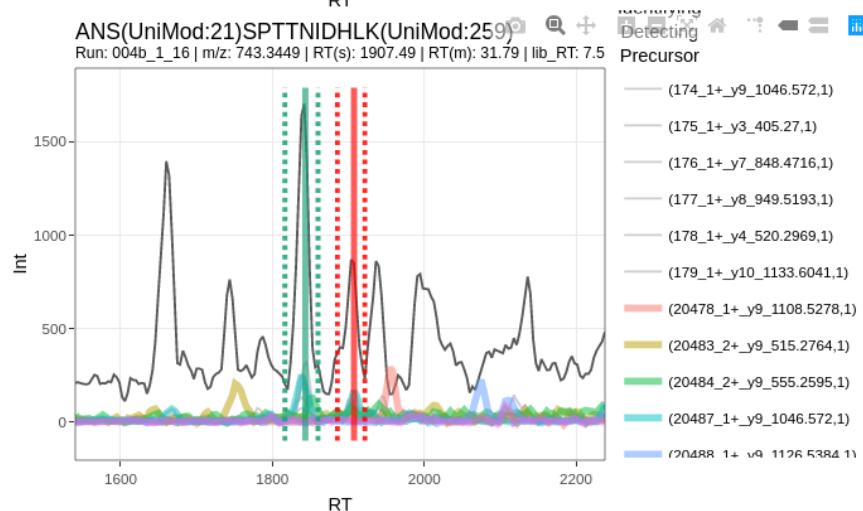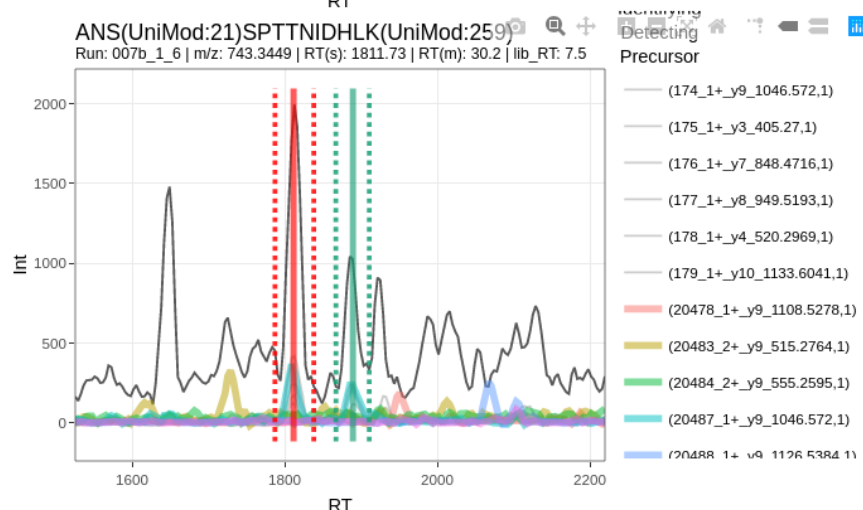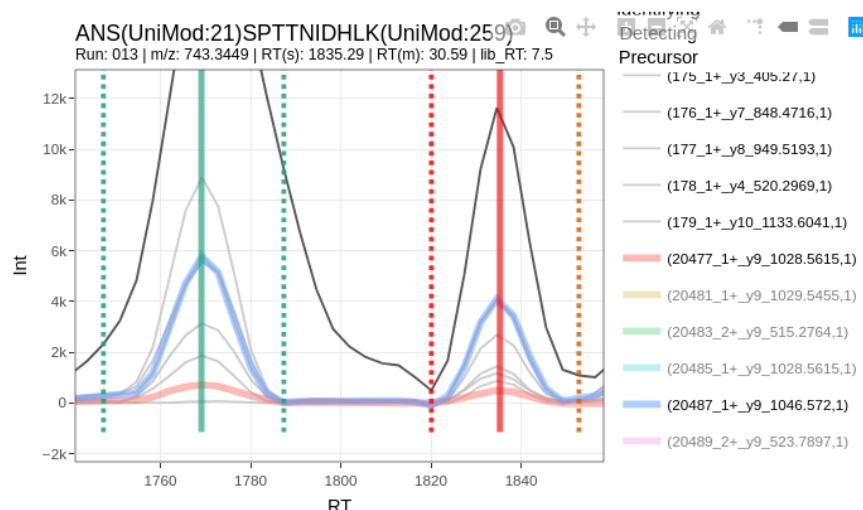

## Displaying a Single Transition

The user can choose to display a single transition by double clicking on the transition legend they wish to display

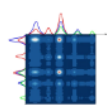

DrawAlignR Ver: 0.1.0

General Settings
Alignment Settings

Plot Settings

☒ Plot Precursor Trace  
☒ Plot Detecting Traces  
☒ Plot Unique Identifying Traces  
Show n Identifying Traces  
  
identifying y-ions  
  
identifying b-ions  
  
☒ Show Transition Scores (hover tooltip)  
☒ Show All Peak-Groups

